



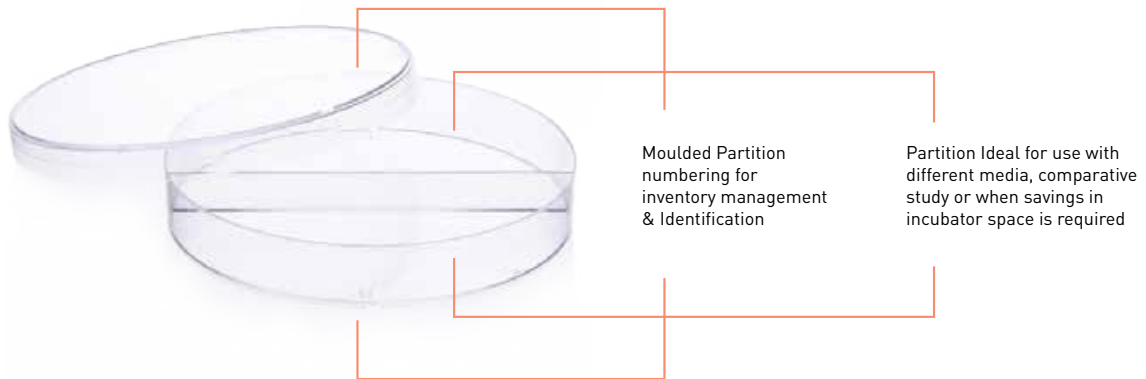

**MICRO BIOLOGY**

2 Compartment  
3 Compartment  
Soft Inoculation Loop  
Cell Spreaders

## PETRI DISH FEATURES

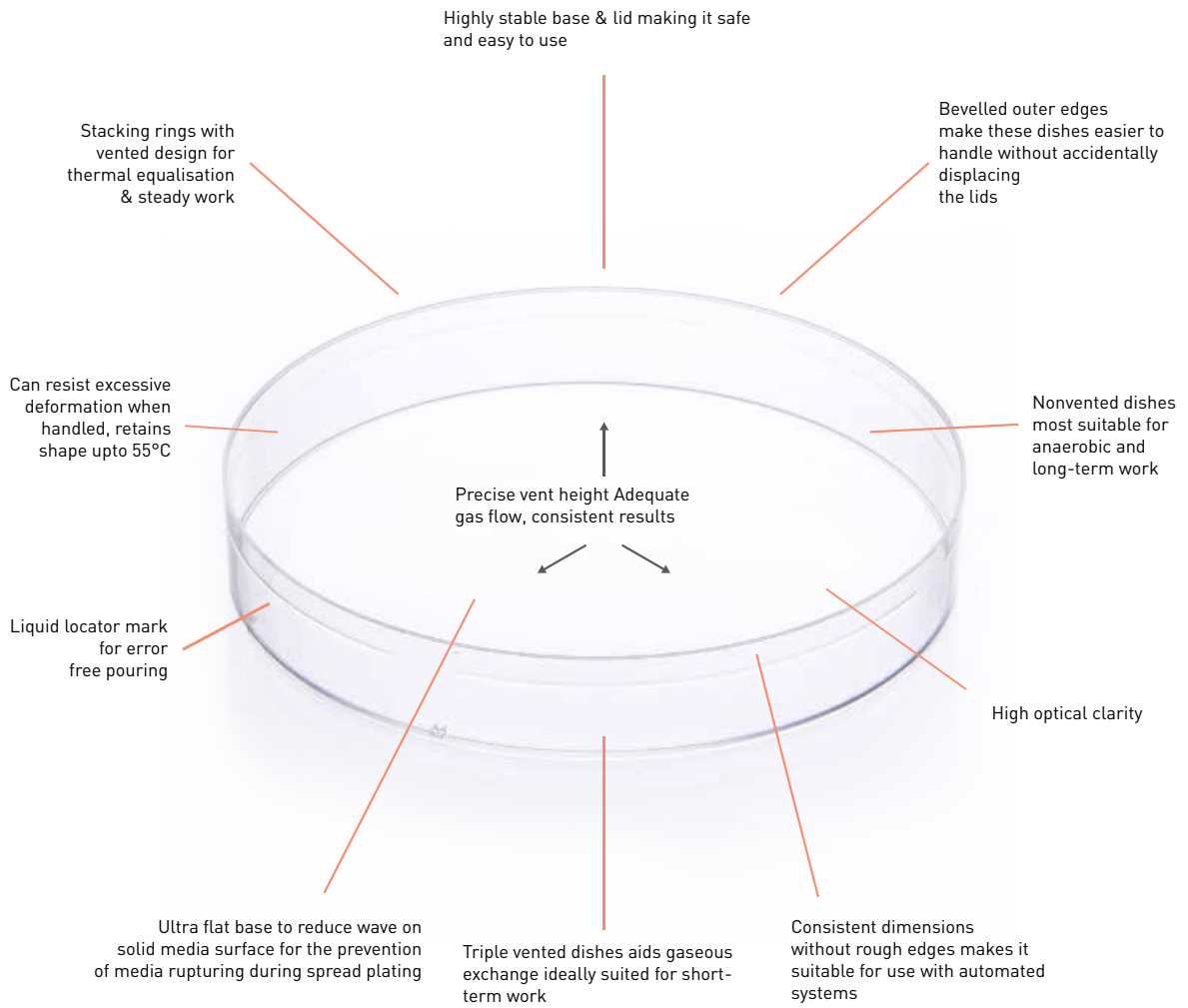

## Petri Dish, PS

### Petri Dish - Aseptic

Code	Diameter (mm)	Height (mm)	Compartment	Vented	Packing	No. per Case	Price / Case (₹)	Price /20 Case (₹)
P10901	90	15.8	No	3 Vent	Sleeve of 20	480	2991	56231
P10902	90	15.8	No	3 Vent	Individual Wrapped	380	2677	50325
P10903	90	15.8	No	0 Vent	Sleeve of 20	480	2991	56231
P10904	90	15.8	No	0 Vent	Individual Wrapped	380	2677	50328

### Petri Dish - Gama Radiated Sterile

Code	Diameter (mm)	Height (mm)	Compartment	Vented	Packing	No. per Case	Price / Case (₹)	Price /20 Case (₹)
P10905	90	15.8	No	3 Vent	Sleeve of 20	480	3312	62258
P10906	90	15.8	No	3 Vent	Sleeve of 10	480	3396	63844
P10907	90	15.8	No	3 Vent	Individual Wrapped	380	3011	56608
P10908	90	15.8	No	0 Vent	Sleeve of 20	480	3312	62258
P10909	90	15.8	No	0 Vent	Sleeve of 10	480	3396	63844
P10910	90	15.8	No	0 Vent	Individual Wrapped	380	3011	56608

### Petri Dish - Gama Radiated Sterile, Double Wrapped

Code	Diameter (mm)	Height (mm)	Compartment	Vented	Packing	No. per Case	Price / Case (₹)
P10928	90	15.8	No	0 Vent	Sleeve of 18	432	3477
P10927	90	15.8	No	0 Vent	Sleeve of 8	384	3262

Abdos high quality disposable Petri Dishes are suitable for all microbiological applications manufactured as per ISO 24998 : 2008 'Plastics laboratory ware — Single-use

Robust dimensional accuracy of 90 mm ensures compatibility with most automated systems.

Extremely high optical clarity for easy view on microscope.

Abdos Petri Dishes retain their shape upto 55° C.

Vented Dishes have elevated lid that allows for superior gas exchange and evaporation.

Non-Vented Dishes have flat lid that prevents evaporation encouraging long incubation period and anaerobic application.

Aseptic Dishes are manufactured in clean room under controlled low bio-burden conditions for best performance in all applications.

Easy to stack.

Shallow dish design allow for maximum use of incubator space.

Abdos Petri Dishes are also available in two & three compartment that enable parallel testing of single sample using different agar, solvents and reagents.

Sterile Plates has been validated according to EN ISO 11137-2; 2012 for SAL 10<sup>-6</sup>



**2 Compartment Petri Dish - Aseptic**

Code	Diameter (mm)	Height (mm)	Compartment	Vented	Packing	No. per Case	Price /Case (₹)
P10911	90	15.80	2	3 Vent	Sleeve of 20	480	3881
P10912	90	15.80	2	3 Vent	Individual Wrapped	380	3465
P10913	90	15.80	2	0 Vent	Sleeve of 20	480	3881
P10914	90	15.80	2	0 Vent	Individual Wrapped	380	3465

**2 Compartment Petri Dish - Gama Radiated Sterile**

Code	Diameter (mm)	Height (mm)	Compartment	Vented	Packing	No. per Case	Price /Case (₹)
P10915	90	15.80	2	3 Vent	Sleeve of 20	480	4302
P10916	90	15.80	2	3 Vent	Individual Wrapped	380	3891
P10917	90	15.80	2	0 Vent	Sleeve of 20	480	4302
P10918	90	15.80	2	0 Vent	Individual Wrapped	380	3891



**3 Compartment Petri Dish - Aseptic**

Code	Diameter(mm)	Height (mm)	Compartment	Vented	Packing	No. per Case	Price /Case (₹)
P10919	90	15.80	3	3 Vent	Sleeve of 20	480	4929
P10920	90	15.80	3	3 Vent	Individual Wrapped	380	4234
P10921	90	15.80	3	0 Vent	Sleeve of 20	480	4929
P10922	90	15.80	3	0 Vent	Individual Wrapped	380	4234

**3 Compartment Petri Dish - Gama Radiated Sterile, PS**

Code	Diameter (mm)	Height (mm)	Compartment	Vented	Packing	No. per Case	Price /Case (₹)
P10923	90	15.80	3	3 Vent	Sleeve of 20	480	5459
P10924	90	15.80	3	3 Vent	Individual Wrapped	380	4919
P10925	90	15.80	3	0 Vent	Sleeve of 20	480	5459
P10926	90	15.80	3	0 Vent	Individual Wrapped	380	4919



## Soft Inoculation Loop Gama Radiated Sterile, ABS

118

Code	Capacity (µl)	Colour	Length (mm)	Packing	No. Per Case	Price / Case (₹)
P10985	1 µl loop + Incorporated Needle	White	200	Sleeve of 20	2000	8500
P10986	1 µl loop + Incorporated Needle	White	200	Individually Wrapped	500	3055
P10987	10 µl loop + Incorporated Needle	Blue	200	Sleeve of 20	2000	8500
P10988	10 µl loop + Incorporated Needle	Blue	200	Individually Wrapped	500	3055

Smooth surface, provide uniform and gentle streaking without damaging the gel surface.


Available in individual wrapped pack and zip lock sleeve that eliminates risk of contamination.

Polygonal shaft improves grip, assists orientation and makes operation easy.


Sterile loops has been validated according to EN ISO 11137-2; 2012 for SAL 10<sup>-6</sup>



STERILE R

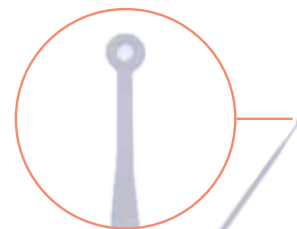

MICRO BIOLOGY

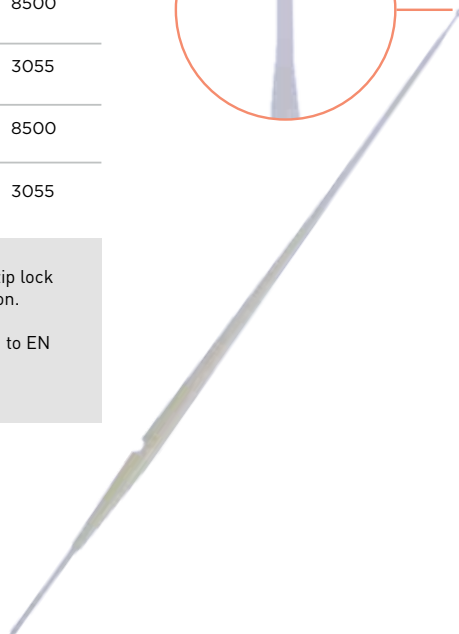

## Cell Spreaders Gama Radiated Sterile, PP

119

Code	Shape	Length (mm)	Blade Width (mm)	Colour	Packing	No. Per Case	Price / Case (₹)
P10989	Triangle Shaped	215	30	Blue	Individually Wrapped	100	1742
P10990	Triangle Shaped	240	60	Yellow	Individually Wrapped	100	2596
P10991	T - Shaped	140	35	Blue	Individually Wrapped	400	5824

Cell spreaders are designed to spread and disperse liquid onto the surface of an agar plate without gouging or cutting the medium.

Sterile Cell Sperader has been validated according to EN ISO 11137-2; 2012 for SAL 10<sup>-6</sup>

Available in Triangular and T-shaped.



STERILE R

