



# STERIKING®

## IMPULSE SEALER IS-250

*Optimal choice for dental & medical clinics and laboratories*

*USER-FRIENDLY, ECONOMICAL*

*ALL-DAY READY FOR USE*

*CONSTANT SEALING PRESSURE, NO WARMING-UP TIME*

*SEALING TECHNIQUE ACCORDING TO DIN 58953, PART 7*

*ELECTRONIC, AUTOMATIC SEALING TIME SYSTEM*

*ACOUSTIC ALARM AND SIGNAL LAMP FOR THE SEALING TIME*

*BUILT-IN CUTTING DEVICE WITH ROTATING CIRCULAR KNIFE*

*PRACTICAL INTEGRATED ROLL HOLDER*

*SAFETY TESTED AND APPROVED BY EN, VDE*

*SHOCK WATCH SYSTEM FOR DELIVERY SAFETY*



**WIPAK MEDICAL**

## IMPULSE SEALER IS-250

With Integrated Roll Holder & Cutting Device

IMPULSE SEALER FOR QUICK AND RELIABLE CLOSING OF HEAT-SEALABLE PACKAGING MATERIALS.


### SIMULTANEOUS SEALING & CUTTING

- Designed for closing all heat-sealable packaging materials in common use.
- Sealing time easily adjustable to suit the material to be sealed.
- Acoustic alarm and signal lamp indicate the adequate sealing time signal can be switched off.
- Fixed constant sealing pressure.



### BUILT-IN CUTTING DEVICE


- Rotating circular knife enables easy manual cutting process in both directions.

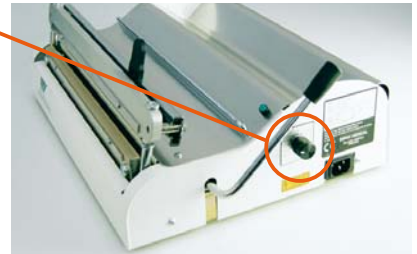

### SEALING TIME SETTING

- Optimal conditions for each material can be set with the adjuster.
- Seals also Tyvek® packages well.
- Not recommended for gusseted packaging.

The IS-250 unit provides a strong seal. Total seal width is 8 mm.

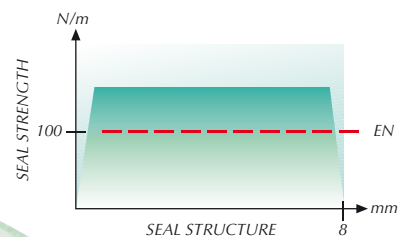

Sealing width max. 250 mm



### TECHNICAL DATA

Sealing time	adjustable
Seal width	8 mm, flat seal
Power consumption	870 W
Power supply	230/115 V, 50/60 Hz
Weight	10 kg
Dimensions	220 x 420 x 300 mm

### WARRANTY

- 24 months from the date of purchase.

### SHOCK WATCH LABEL

- Enables recognition of transport damage.
- Adhered to the transport carton.

Detailed Operating Manual and List of Spare Parts come with each unit.



**WIPAK MEDICAL**

Wihuri Oy Wipak  
P.O. Box 45, Wipaktie 2  
FI-15561 Nastola, Finland  
Tel. +358 (0)20 510 311  
Fax +358 (0)20 510 3333  
E-Mail: info.finland@wipak.com

www.steriking.fi www.wipak.com